




 EFFICIENT IS **SILGINAT**[®]

ANATOMIC IMPRESSION

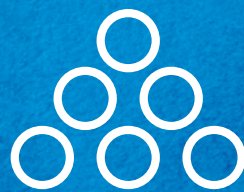

KETTENBACHDENTAL

Simply intelligent

SILGINAT® – SIMPLE, FAST, EFFICIENT!

Anatomical impression

Silginat® is an alginate replacement material based on A-silicone especially for **alginate applications**.









INDICATIONS

- **Opposing dentitions**
- **Fabricating temporary crowns and bridges**
- **Anatomical models**
- **Fabricating simple removable prosthetic restorations**
- **Producing models for the construction of trays for whitening solutions, bites and surgical guides**
- **Orthodontic work**
- **Case study models**
- **Fabricating removable retainers and splints**




Technical data

 Working time at 23 °C / 74 °F ≤	 Intraoral setting time at 35°C / 95 °F ≥	 Total setting time* ≥	 Elastic recovery test (approx.) %	 Linear dimensional change (approx.) %	 Hardness (ap- prox.) Shore
90 sec.	90 sec.	3:00 min.	99.5	– 0.20	A 45

*Total setting time (removal time from mouth) from the start of mixing

AN A-SILICONE AS YOU KNOW IT

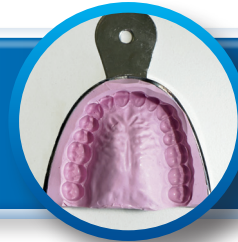

- Storage stability and dimensional stability, high precision

Silginat vs. alginate
approx. 24 h after impression-taking



EASY TO HANDLE



- Clean and fast work thanks to the 5:1 jumbo cartridge with reproducible results



COST-EFFECTIVE IMPRESSION TAKING

- Pourable multiple times, only one anatomical impression per patient required



FOR EVEN GREATER EFFICIENCY
Saves time and many work steps during your work!



SILGINAT® – THE RANGE

Silginat®		
	Silginat®	Silginat® Strawberry
Refill pack	2 x 380 mL cartridges	2 x 380 mL cartridges
Normal pack	6 x 50 mL cartridges, 6 mixing tips	6 x 50 mL cartridges, 6 mixing tips
Bulk pack	24 x 50 mL cartridges	24 x 50 mL cartridges



NOTES

Kettenbach GmbH & Co. KG
 Im Heerfeld 7
 35713 Eschenburg, Germany
 Phone +49 (0)2774 7050
 info@kettenbach.com
 www.kettenbach-dental.com


KETTENBACHDENTAL
 Simply intelligent



027745_4523