





New Cartridge System

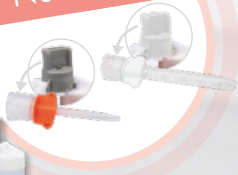

# SOOO ADHERENT IS MUCOPREN® SOFT

SOFT RELINING

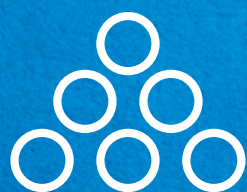


KETTENBACHDENTAL  
Simply intelligent




# MUCOPREN® SOFT

## Prosthetic relining

Mucopren® Soft is a soft long-term relining material based on a silicone for **direct and indirect application**.



**Special long-term adhesion:**  
Excellent adhesion to acrylic dentures through patented Mucopren® Adhesive



**Unique sealing through its very smooth, hydrophobic silicone surface:**  
Inhibits bacterial colonization thus ensuring hygienic conditions on the prosthesis<sup>1)</sup>



Step by Step

### Simple treatment & fast processing

Faster setting times and shorter waiting times ensure a fast and smooth process in the practice and laboratory.

### Supports the healing process of implants

Also suitable for use in implantology for trouble-free osseointegration.



### High wearing comfort

The definitive resilient reline is completely odorless and tasteless and makes wearing the prosthesis as pleasant as possible.

# HOLDS TO ITS PROMISES



## FURTHER INDICATIONS



- For the creation of **permanently soft relinings**
- For relieving pressure on the surgical site after insertion of **implants**
- For balancing mucous membrane areas with different resilience
- For relieving sharp bony ridges
- For relining **interim prostheses**
- As **mucous membrane conditioner** in case of prosthetic stomatitis
- For creating and relining **obturators**

## Benefits of the optimized adhesive<sup>2)</sup>

» ONLY A SINGLE  
APPLICATION  
NECESSARY



3 MIN

### CHAIRSIDE

The processing time is reduced to only **3 minutes** following mouth removal. No water bath necessary.



3 MIN

### LABSIDE

The processing time is reduced to only **3 minutes** in the water bath or pressure pot.

<sup>1)</sup> Result from in-vitro studies in which it was shown that the growth of *Candida albicans* was prevented thanks to the smooth silicone surface

<sup>2)</sup> Refers to the optimized product version of Mucopren® Adhesive, first introduced in 2023.



# MUCOPREN® SOFT – THE RANGE



## Mucopren® Soft

<b>Mucopren® Soft Basic Set</b>	REF 28105	50 mL Mucopren® Soft, 50 mL Mucopren® silicone sealant, 10 mL Mucopren® Adhesive, 6 mixing tips white, 16 mixing tips orange, 1 brush holder, 20 single-use brushes, 1 steel bur, accessories
<b>Mucopren® Soft Normal pack</b>	REF 15687	2 x 50 mL Mucopren® Soft, 6 mixing tips white
<b>Mucopren® silicone sealant</b>	REF 15686	50 mL Mucopren® silicone sealant, 8 mixing tips orange
<b>Mucopren® Adhesive</b>	REF 14203	10 mL Mucopren® Adhesive



### Technical data

Product name	Mixing ratio	Working time at 23 °C / 74 °F ≤	Intraoral setting time at 35 °C / 95 °F ≥	Water bath	End of setting ≥	Shore hardness (approx.)
Mucopren® Soft – direct relining	1:1	2:15 min.	3:15 min.	Not applicable	5:30 min.*	A 28
Mucopren® Soft – indirect relining	1:1	2:15 min.	Not applicable	3:00 min.	> 3 min.	A 28
Mucopren® silicone sealant	1:1	1:30 min.	Not applicable	Not applicable	6:30 min. at 23 °C / 74 °F	A 33

\*Total setting time (removal time from mouth) from start of mixing

## NOTES

---



---



---



---



---

Kettenbach GmbH & Co. KG  
 Im Heerfeld 7  
 35713 Eschenburg · Germany  
 Telefon +49 (0)2774 7050  
 info@kettenbach.com  
 www.kettenbach-dental.com



1410871\_4924